

8 in 1 Multifunction skin care machine

SPA10

User Manual



Warning

Dear user:

Thanks for using our company's products, in order to fully make use of the products, we sincerely suggests follow:

◆◆Please read the instruction in detail and keep it well for further reading and consulting.

◆◆Please follow the instruction to fit and operate the machine correctly.


◆◆Please do not remove and change any accessories of the machine.

◆◆It forbids others to open and remove the mainframe except the authorized people in the company.

◆◆It forbids others to change and maintain the machine except the professional people.

◆◆Please contact us if there is any trouble during operation.
This product does not belong to medical, please don't use for disease treatment!


The back of the machine:

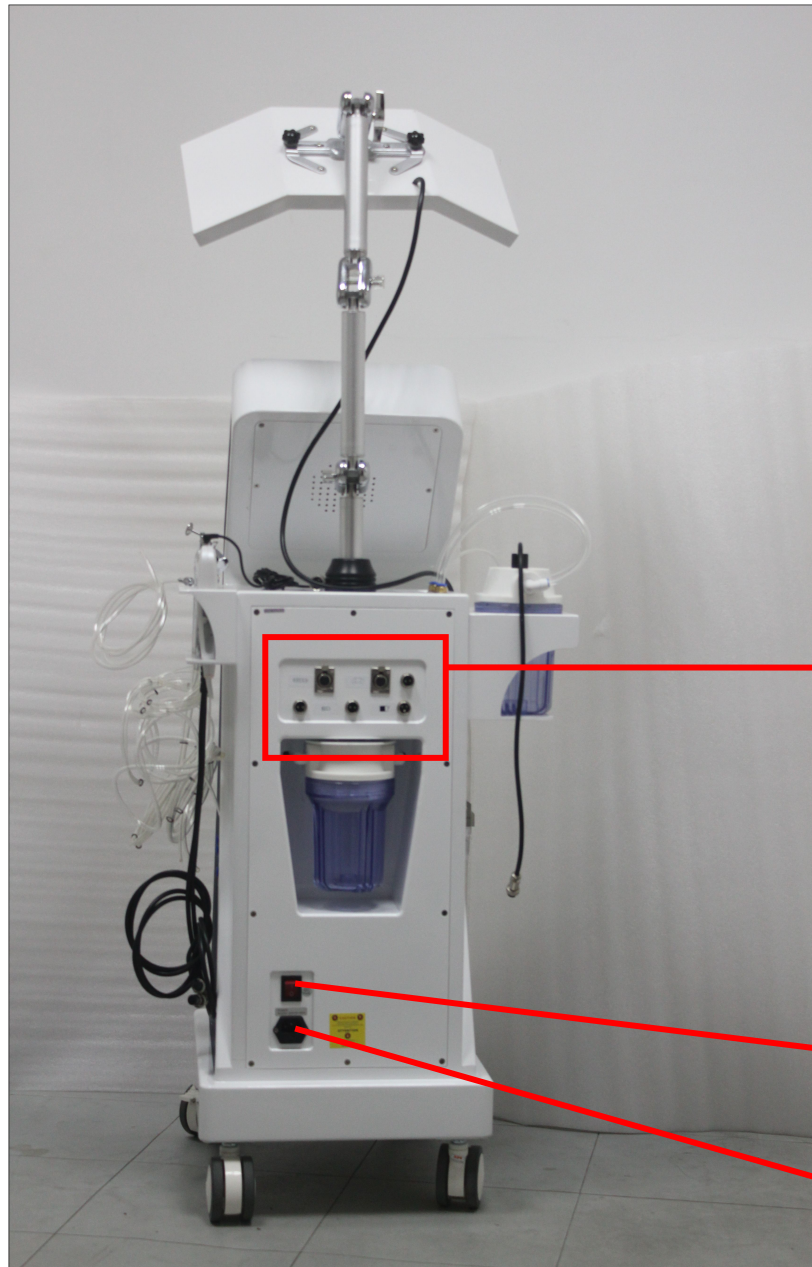

PDT LED

skin scrubber

Recycle water bottle

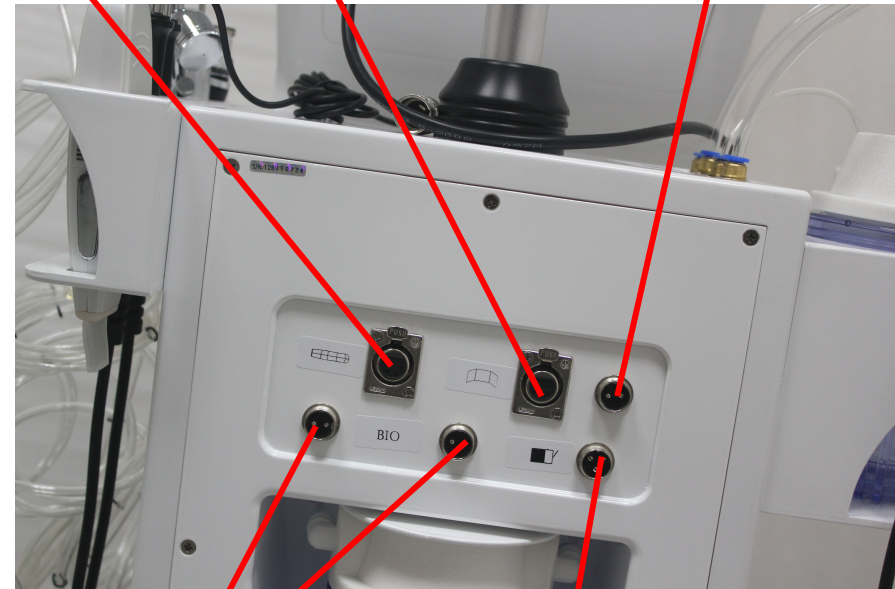


BIO

Foot Switch

General power switch

Power supply line jack

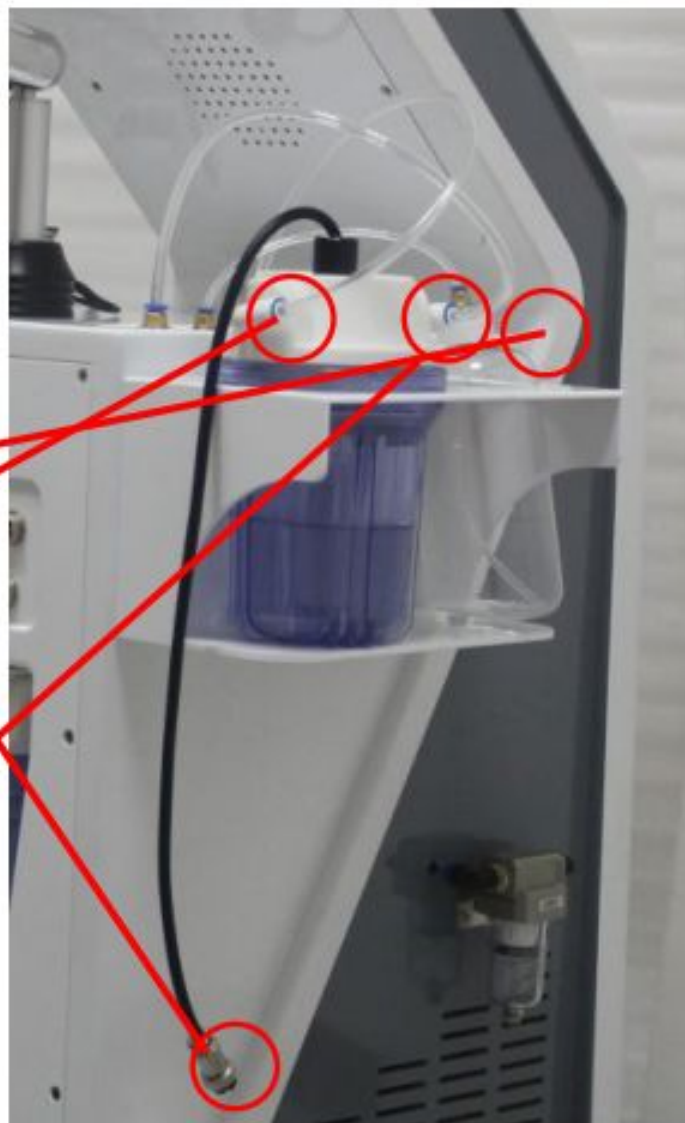

Parts Introduction:



Hydro dermabrasion



Diamond dermabrasion



Oxygen sprayer



Oxygen injection



PDT LED



Skin scrubber



BIO



Water oxygen

Basic settings:

Step 1: Connect the accessories to the host.


Step 2: Insert power line, open the General power switch

Step 3: Input the default password "000000" picture<1>

Step 4: Enter interface as. Click "basic setting" to setting language. picture<2>

Notes "advanced setting" is for manufacturer only.


Step 5: Choose the language, we have English and Chinese. picture<3>

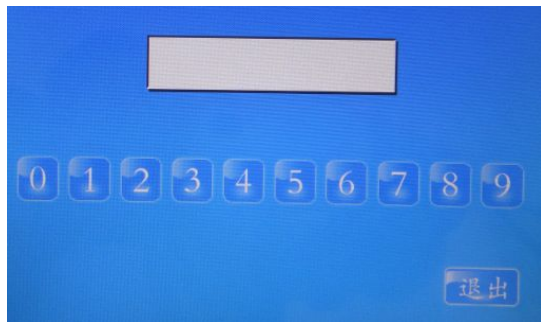

Picture<1>



picture<2>

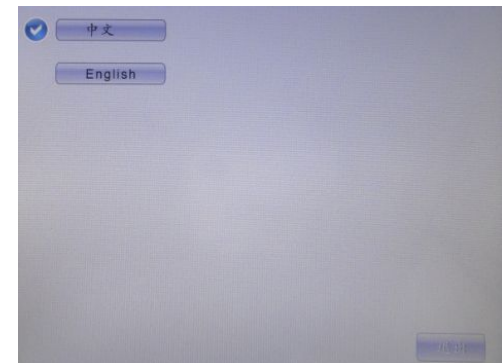

picture<3>

Principle of hydro dermabrasion:

Hydro dermabrasion is an exciting NEW technology that combines the efficacy of Micro dermabrasion, vacuum system and new Aqua Fuse hydration system.

It is a lot gentler than crystal micro dermabrasion or diamond dermabrasion as only water is used, and we use professional quality machines.

Hydro dermabrasion (or hydra micro dermabrasion) applies mechanical and chemical peeling simultaneously.

A hydro-dermabrasion machine includes a water jet, an air compressor, a two-way flow control valve, a purified container and a waste water tank.

Unique hydra-dermabrasion tip jets out thin and tiny stream of water and polish skin in high speed by vacuum suction, the liquid from which helps moisture skins in the mean time.

Hydro dermabrasion operating steps:

Step 1: Power line and aqua peel handle connect to mainframe.

Step 2: Vortex tips connect to aqua peel handle.

Step 3: Put essence into clean water bottle and then connect clean and recycle water bottle to mainframe.

Step 4: Turn on the power and put handle onto face skin.

Step 5: Move evenly as picture and adjust water strength properly.

Step 6: 10-15 minutes operation time control is better.

Step 7: You should take pipe out correctly and brush tips after operation.

Principle of diamond dermabrasion:

Diamond dermabrasion operates on the premise that encouraging exfoliation of the skin is desirable.

The word “diamond” is used to separate this procedure from true dermabrasion, where the deeper layers of the skin are removed.

By using a motorized brush system or bombarding the skin with high-speed particles, it is possible to remove desquamating corneocytes, and make the skin smoother and softer with a better color.

It is important not to go any deeper into the skin than the first few layers of stratum corneum to avoid scarring and other problems.

Diamond dermabrasion can be used by trained skin care professionals because it does not enter the living layers of the skin, any procedure that enters the living layers of the skin and causes bleeding should be performed by a physician.

Diamond dermabrasion operating steps:

Step 1: link power cord and glass metal handle to mainframe.

Step 2: put water into clean and connect recycle water bottle and clean water bottle to the machine.

Step 3: pipe and diamond tips connect to handle and press down the power switch.

Step 4: put the glass metal handle tightly onto corresponding body place according to different tips and you can see water flowing via pipe.

Step 5: move evenly as picture and adjust strength of water.

Step 6: the best operation time is 10-15 minutes one time.

Step 7: you should take pipe out correctly and brush tips after operation.

Effect:

Hydra facial dermabrasion help to improve the appearance of

- Fine lines & wrinkles
- Hyper pigmentation & melasma
- Sun damaged skin
- Age spots
- Blackheads & whiteheads
- Enlarged pores
- Excess oil
- Hydration
- Puffy eyes
- Rejuvenation,over all skin health
- Removes skin debris, scars, blemishes
- Removes wrinkles and uneven pigmentation of the skin

Notice:

Don't put any liquid onto or into the machine. Be careful about the electric connections.Avoid operations on pregnant women or those people with epilepsy, pacemakers or any serious psychological or medical condition.Please disassemble or assemble the equipment by professional personnel.

Principle of Oxygen spray & injection:

Oxygen is indispensable to the human body's "air of health," is the energy needed for body as a source, it is the first needs of human life. oxygen therapy, developing rapidly, both in the development history of medical or cosmetic, oxygen therapy has made a significant contribution, very well received, has generally been recognized as a widely used technology of medical beauty.

Almighty Oxygen instrument uses the world's advanced PSA technology (PSA), at room temperature, after connecting the electricity, it can isolate from the air, and has the effect of atomized jet of oxygen, injecting Oxygen, inhaling oxygen, solve the problem that short of oxygen of. important the continuity of high concentrations of dust-free, sterile medical oxygen through body, skin fully Almighty Oxygen instrument applies foreign advanced technology, non-needle injection of oxygen through a variety of specially designed oxygen output tools push the pure oxygen and products containing active elements into deep skin through natural gap between the epidermal cells make skin quickly absorb oxygen, enhance activity of human cells and speed up cell metabolism, so that skin regains its healthy, moist, delicate young state be charming and glory.

Oxygen sprayer operating steps:

- Step 1:** Link power cord and spray handle to mainframe.
- Step 2:** Connect pipe to spray and turn on the machine.
- Step 3:** Add water to spray.
- Step 4:** Start spray and you can see the water mist ejecting.
- Step 5:** Put spray handle onto face and then move evenly on the face.
- Step 6:** You can adjust the pressure intensity according to your adaptive capability.
- Step 7:** Operation time is better 3-5 minutes doing like this.

Oxygen injection operating steps:

- Step 1:** Link power cord and spray handle to mainframe.
- Step 2:** Connect pipe to injection and turn on the machine.
- Step 3:** Put spray handle onto face and then move slowly on the face.
- Step 4:** Operate for 7-10 minutes.

Effect:

oxygen sprayer together with the active concentrated essence, to do the pure oxygen massage on the skin in order to treat all kinds of prurituses, sensitive skin, neurotic dermatitis, contact dermatitis, eczema, wheek, the sunlight symptom, etc.

Helps to calm down the skin and to enhance the healing processes of the tissues, it can be used also in case of active acne, acne rosacea and herpes simplex, in order to calm down the inflammatory process
oxygen injector treat and improve the problems that caused by the oxygen shortage, such as dark and light skin, wrinkles, looseness, pigmentation and the conventional skin problems.

Principle of PDT LED:

Beam Therapy was brought forward early in 18 century, at the same time , the recipe to cure tuberculosis is of sunshine shining.

In 1903, a professor in Denmark was awarded Nobel Prize for his research of utilizing ultraviolet to cure phthisic skin as well as red light to cure small pox.

The instrument makes use of power light caused from 624nm pulse to procure very excellent therapeutic effect without any side effect to eyes. It doesn't hurt or break skin. Power light shining scarfskin can stop decomposition and cohesion of black pigment.

When power light is shining healthy skin, the energy will decrease automatically, and collagen in hypodermal organization will realign, with view to enhance the luster and spring, contract skin pore and lighten small wrinkle.


PDT LED operation steps:

Step 1:Put all accessories well and connect to the host.


Step 2:apply the gel on the face then remove in circle.And wear the glasses

Step 3:Set the "PDT LED model"

Step 4:you can set working time firstly,in general,the time is 30 minutes,according to your requirment set working time(reduce or increase working time)

Step 5:set"  "you can choose different colors.

Step 6:set"  "you can choose differnt intensity.

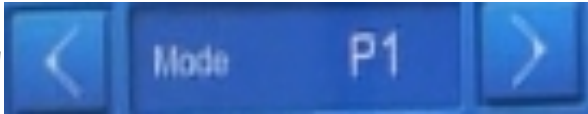Step 7:set"  "you can choose different model,it has different frequency with LED light.

Step 8:Times out,the machine automtic.

Effect:

Red light:Promote blood circulation, activate muscles, reduce fat, relieve anemia and shortness of breath, improve cardio-pulmonary function, and evenly largely benefit mood.

Blue light:Promote the synthesis of protein and gastric collagen, activate muscles,contract the slack skin and reduce lines of pregnancy. Irradiation to the point of glabella can compose mood and relax body and soul.

Yellow light:Enhance blood circulation,activating cells, stimulate the cell excitability active role. Increasing the skin collagen fiber and elastic fibers, avoiding the skin redness.

Green light:Has the effects of neutralization,balancing and calming,it can ease strain,dredge lymph and useful in edema elimination.

Red+Blue PDT:cure the acne, relieve skin.

Red+Yellow PDT:Stimulate dermal collagen, improve skin lines and skin tone .

Green +Blue PDT:Improve facial lymph circulation, inhibit pigment expansion .

Red+Blue+Green PDT:Improve ability of the cells absorb, promote collagen formation and treat the skin.

Principle of skin scrubber:

With Ultrasonic Scrubber begins a new generation of aesthetic soft peeling.

Using the high performance of ultrasounds, this exciting cosmetic breakthrough provides a safe painless and perfect physiological skin regeneration.

By removing the superficial scaly layers of the skin, the Scrubber stimulates the revival of the skin tissue rendering it youthful and luminous from the first treatment.

Skin scrubber operation steps:

Step 1: Link power line to mainframe.

Step 2: Connect skin scrubber to mainframe.

Step 3: Turn on power.

Step 4: Put scrubber onto skin and Must put your finger on both sides of skin scrubber as picture (skin scrubber it has 4models can change adjust on the handle).

Step 5: Move slowly and evenly on the skin as follows.

Step 6: The best operation time is 10-15 minutes.

Effect:

- 1.Deep clear and clean up aged horniness.
- 2.Reduce melanin to get a whiten skin.
- 3.Deep fine wrinkle and make skin young and soft.
- 4.Improve blood circulation and lymph metabolism.
- 5.Make skin immediate glow and regrow.
- 6.Accelerate the absorption of nutrition.
- 7.Digital control.
- 8.No side effect for any type of skin.
- 9.Slick design and convenient operation, easy and small enough to keep.

Principle of BIO:

Skin is the body's largest organ, composed by the multi-layer structure, while the collagen layer which decides skin softness and flexibility will become aging as people get old.

BIO wrinkle removal simulates the biological current issued by the brain, supported with a professional operation technique to lift up skin with different levels speed up the skin's blood circulation, increase metabolism and promote skin cell regeneration and repair of collagen to help restore the skin's natural elasticity, effectively remove wrinkles, eliminate dark circles.

BIO operation steps:

Step 1: Deep clean the skin.

Step 2: Apply some wrinkle removal products (essence, gel) on skin need cared.

Step 3: Set work time: Face (including neck) care about 30 minutes, eye-care about 20 minutes.

Step 4: Adjust current output intensity (1 ~ 60), current output should be from weak gradually to strong, within acceptable range of customers.

Step 5: Beautician holds 2 wrinkle removal handle, each hand hold one, to do wrinkles removal, skin lifting care, the two wrinkle removal hand should touch skin, but do not stay on any part over 8 seconds.

Step 6: Clean facial skin with warm water. means location point for wrinkle removal handle means move direction

Principle of Water oxygen:

Use of hyperbaric oxygen, push the sterile saline (saline) and oxygen access air compressors simultaneously, transfer water drops of the particle type into small molecular (5-200 μ m microns), tiny water-drops would spray (jet spray), with 200m/sec (seconds) speed (supersonic) sprayed by a specially designed nozzle, inject water and oxygen into the deep skin speedy.

Exfoliating and do needle-free nutritional delivery treatment painless noninvasively.

The impact generated by hyperbaric oxygen and saline, physical exfoliating procedure suit for any skin, it can stimulate the cells regeneration function, stimulating blood circulation, strengthen the metabolism of skin cells and collagen proliferation!

Notice: use the saline water

Water oxygen operation steps:

Step 1: Link power line to mainframe.

Step 2: Deep clean the skin.

Step 3: Adjustment the intensity suitable with you.

Step 4: Holding water oxygen spray pen, stepping on foot switch, then spray at the face.

Step 5: Stop to step foot switch, water oxygen end.

Effect:

1. Removal of scar: burn surgery scar acne pits, Pockmark and so on.
2. Removal of acne: remove acne and inflammation, acne, blackheads, acne and acne scar.
3. Skin rejuvenation: remove splash, Pigment deposition, shrink pores, whiten skin, enhance skin elasticity and luster improving skin all over the body.
4. Wrinkle removal: Lips wrinkles, eye wrinkles improve glabellar lines, stretch marks and so on.
5. Hair growth: improve seborrheic alopecia alopecia areata, hair care.
6. Accelerate blood circulation: increase skin elasticity, reduce skin relaxation.
7. Repair damaged skin.

Safety notice:

- 1.The machine is professional instrument for beauty salon, in order to avoid any unnecessary hurt, need professional beautician to operate
- 2.Do not disassemble the machine or operate the machine that not explained in the instruction book. All the maintenance services shall be carried out by the professional maintenance technician approved by our company .
- 3.Do not install or operate the machine with wet hands or near the water sources so as to avoid the liquid spraying/ swashing over the machine.
- 4.Pull out the plug and connect with the maintenance service center in the following cases:
 - A.liquid spray/swash over the machine or he machine get caught in rain or water .
 - B.the machine has special smell, fume/smog or abnormal noise.
 - C.power lines damaged or scratched.
 - D.instrument fell down or the shell damaged.
- 5.Do not put any item on the power wires, or install the instrument in the place where its power wires will be trampled.
- 6.Pull out the power wires at the end of use for the sake of safety.
- 7.Do not insert any type of item into the instrument opening; otherwise there will be a fire or electric shock. In case of something fell into the opening but you could not take it out, please pull out the power wires and connect with the dealer or our company.

Technique parameter:

Working frequency:50Hz

Electricity voltage:AC 11V/220V 50-60Hz

Power:800W

Wave length of PDT:Red Light $630\pm 5\text{nm}$

Blue Light $460\pm 5\text{nm}$

Treatment Intensity of PDT: 250mW/cm^2

Vacuum:80 Kpa

Accessories List:

| | |
|------------------------|------|
| 1.Host | 1pc |
| 2.Annex shelf | 1pc |
| 3.Power supply line | 1pc |
| 4.Foot switch | 1pc |
| 5.Water dermabrasion | 8pcs |
| 6.Diamond dermabrasion | 9pcs |
| 7.Oxygen sprayer | 1pc |
| 8.Oxygen inject | 1pc |
| 9.PDT LED | 1pc |
| 10.BIO | 2pcs |
| 11.Skin scrubber | 1pc |
| 12.Water oxygen | 1pc |
| 13.User maunal | 1pc |